BLACKROLL®

BOOSTER

Instruction Manual

EN | DE | FR | NL | IT | PL | CN


Safety Instructions

- + Do not use on sensitive areas or near organs.
- + In case of pain or other contraindications, stop use immediately and consult a doctor.
- + Persons with reduced physical, sensory or mental capabilities and any person under the age of 18 must be supervised during use.
- + Do not insert electrically conductive objects into the openings of the appliance.
- + Do not insert the appliance into any orifices.
- + Never block the vibration of the device. The device must be able to vibrate freely.
- + Do not use as a support or substitute for medical applications without consulting a doctor.
- + For all applications in the area of the head, neck and upper back, use the vibration carefully and selectively.
- + Do not immerse the device in water.
- + Keep liquids away from the ventilation slots and the control buttons.
- + Keep the unit away from heat sources, children and pets.
- + Protect the unit from moisture, dust, aggressive liquids and vapours.
- + If moisture enters the unit, switch off the unit immediately and, if necessary, disconnect the mains plug from the socket.
- + Do not remove any screws or attempt to disassemble the unit.
- + Do not pass the unit on to third parties without the operating instructions.
- + In case of misuse, all warranty claims are void.

Contraindications

Basically, there are hard indications that exclude the use of the device as well as contraindications for which an individual decision must be made as to whether the device can or may be used. If one or more of these points apply to you, please consult your doctor first:

- + Redness
- + Swelling
- + Acute pain
- + Osteoporosis
- + Damage to intervertebral discs
- + Taking anticoagulant (blood thinning) medication
- + Open, inflamed or infected skin surfaces
- + Fresh fractures
- + Increased bleeding tendency, tendency to haematoma formation (formation of bruises)
- + Tendency to embolism
- + Cardiac pacemakers (HSM, PM)
- + Direct bone or eye contact
- + Fibromyalgia
- + Rheumatic diseases
- + Joint replacement
- + Tumour diseases
- + Pregnancy
- + Diabetes

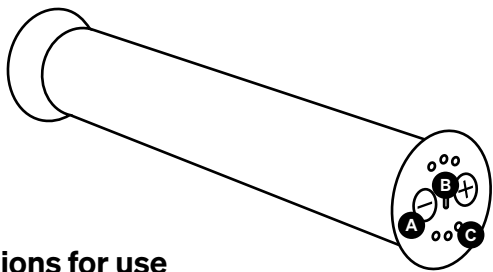Instructions for use

- A** On/Off button & change of vibration frequency
- B** Micro-USB connector
- C** Vibration frequency & battery status indicator*

1. Check all parts for damage.
2. Fully charge the battery using the micro-USB connector and the supplied cable before first use. For safety reasons, do not use the unit while charging. Charging time: approx. 7-10h. Please use the 5V 2A charger when charging, please DO NOT use any fast charger.
3. To switch the unit on or off, press and hold the + and - buttons simultaneously for 3 seconds.
4. When switched on, you can change the vibration frequency by briefly pressing the +- buttons.

The device switches off automatically after 10 minutes. If you still want the vibration, switch the unit on again by pressing the +- button.


*Permanent light = min. 90% charged. At 100%, the charging process stops automatically; the LED remains lit. Shortly after switching on or by pressing a button during the charging process, the charging status can also be displayed via the 6 LED lights.

Vibration frequency 12-56 Hz

The unit covers a wide and precisely controllable vibration spectrum. The frequency range is 12-56 Hz. The deep-acting vibration can irritate the pain receptors. This can lead to a reduction in the sensation of pain. The blood circulation of the muscles can be stimulated, which can promote regeneration.

When do I use which frequency?

We recommend using the low frequencies up to 30 Hz for regeneration and release of muscle tension, e.g. after training. The frequencies above 30 Hz are particularly suitable for activating the muscles, e.g. before training.

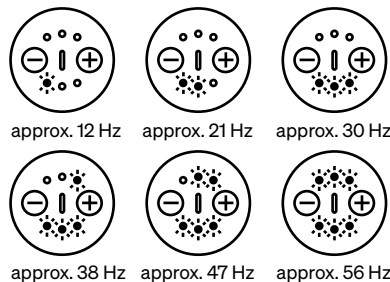Frequency range BOOSTER: 12-56 Hz

The LED lights indicate the intensity of the vibration.

Areas of application

The BLACKROLL® BOOSTER can be used in 3 different ways and thus offers extensive possibilities of use.

1. application in the fascia roller

The BOOSTER can be combined with any 30 cm long BLACKROLL® fascia roller. To do this, unscrew the cover with lettering on the BOOSTER and slide the BOOSTER into a 30 cm long BLACKROLL® fascia roller. Align the BOOSTER and firmly screw on the cover with the lettering on the opposite side. Set the BOOSTER to the desired frequency and start self-massage. The self-massage can be intensified and made more effective by the deep-acting vibrations.

2. use as a stand-alone device

The BOOSTER can be used individually in functional vibration training. By combining exercises with the vibrating bar, the deep muscles can be strengthened effectively.

3. use with BOOSTER HEAD attachments

The BOOSTER can be combined with the BOOSTER HEAD for spot vibration massage. The BOOSTER HEAD includes 4 different attachments in different degrees of hardness and shapes with which pain points can be specifically treated.

Exercises with the BOOSTER.

Download the app now for free and start.

Cleaning

- + Do not clean the unit while it is charging.
- + Do not use water to clean the unit, only a damp cloth.
- + Do not use corrosive cleaning agents or solvents.

Technical specifications

Plastic vibrating core

Vibra Motion technology: eccentric vibrations

Developed in Germany, CE certified

Lithium-ion battery:

Power supply:

Input power supply:

Output power supply:

Max. Power consumption:

Battery:

Charging time:

Battery life:

Charging via micro USB port

6244-FW8002.1/USB

(power supply available on request, not included)

Please use the 5V 2A charger when charging, please

DO NOT use any fast charger.

100-240V / 50-60Hz

5V / 1.4A

45 W

3x S1P = 11.1 V 2000 mAh

approx. 7-10h

approx. 2h at highest level

Disposal

The unit contains electronic components and a rechargeable battery. It must not be disposed of in normal household waste. All electrical and electronic equipment must be disposed of at government designated sites.

Disclaimer

The unit is used for exercise at your own risk. It is recommended that you seek advice from a physiotherapist, manual therapist, masseur or medical specialist if you experience severe pain before or after self-massage or recurrent injuries.

Warranty

The seller grants the buyer, if he is a consumer, the statutory warranty in accordance with the legal requirements for a period of 2 years from the purchase of the item. The Seller's warranty obligations shall initially be limited to supplementary performance, whereby the Seller shall be free to choose whether to deliver a defect-free item or to repair the defective item. If subsequent performance is unsuccessful, the buyer shall be entitled to the further statutory claims. The seller rejects warranty claims for damage caused by accident, obviously incorrect use, misuse or improper transport; which occurred as a result of repairs carried out by unauthorised personnel; which occurred as a result of use and handling not in accordance with the operating instructions; which occurred as a result of improper handling of the purchased item. Furthermore, claims for damages by the purchaser are excluded. Excluded from this are claims for damages by the buyer arising from injury to life, limb or health as well as liability for other damages based on an intentional or grossly negligent breach of duty by the seller, his legal representatives or vicarious agents. Also unaffected is the liability for breaches of duty, the fulfilment of which makes the proper execution of the contract possible in the first place and on the observance of which the buyer may regularly rely. In the event of a slightly negligent breach of such contractual obligations, the Seller shall only be liable for the foreseeable damage typical for the contract, unless it is a matter of claims for damages by the Customer arising from injury to life, limb or health.

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Get the BLACKROLL® App. Free download.

Discover how to get the most out of your BOOSTER. Find exercises and inspiration in our free app.

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Déchetterie

This product complies with the requirements of Council Directive 2014/30/EU on electromagnetic compatibility.

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EN | DE | FR | NL | IT | PL | CN

Product description

Congratulations on your new BLACKROLL® fascia product! These fascia rolls and balls made of EPP (expanded polypropylene) are perfect for self-massage and fascia training. They promote blood circulation, release tension, relieve pain and improve mobility. The products are available in different sizes and surface structures and can be used for different parts of the body as needed.

Sizes and degrees of hardness

Roll length:	6 cm, 15 cm, 30 cm, 45 cm
Roll diameter:	3 cm, 5 cm, 10 cm, 15 cm
Degrees of hardness:	soft (80%), medium (100%), hard (150%)
Ball diameter:	8 cm, 12 cm

Surfaces

Smooth:	all-rounder for regeneration and activation.
Saw blade:	textured surface for increased blood flow.
Grooved:	activating effect with intense increase in blood flow.

Safety instructions

- + Use the fascia products on a non-slip surface to avoid slipping or falling.
- + Pay attention to the controlled application of pressure: Excessive pressure can cause bruising or tissue damage.
- + Consult a doctor if you suffer from health problems such as osteoporosis, blood clotting disorders, herniated discs or varicose veins.
- + Avoid application to acute injuries, redness or inflammation.
- + Store the products in a dry place and not in direct sunlight.

Contraindications

Do not use the BLACKROLL® fascia rolls and balls if you have osteoporosis, blood clotting disorders, fresh injuries or inflamed skin areas without consulting your doctor first. The products should not be used on areas with acute pain or inflammation.

Application

1 Preparation

- + Make sure that the surface on which you are using the fascia products is stable and non-slip.
- + Start with light pressure and increase it if necessary.

2 Use

- + Slowly roll the fascia roll or fascia ball over the affected muscle area.
- + Make sure to adjust the intensity of the massage to your well-being.
- + Smaller rolls and balls are better suited for targeted treatments of small muscle groups (e.g. feet, arms).

3 After use

- + Clean the products with water and a soft cloth.
- + Store them in a dry place to ensure the longevity of the material.

Technical data

Material:	EPP (expanded polypropylene)
Hardness:	Varies by product (soft, medium, hard)
Insoluble in water:	Yes
Recyclable:	Yes

Care instructions

- + Clean the fascia products regularly with water and a soft cloth.
- + Avoid aggressive cleaning agents or chemical substances that could damage the material.
- + Store the products in a dry place and protected from direct sunlight.

Possible uses

- + Self-massage: Use the fascia rolls and balls for targeted massage of muscle groups such as the back, neck, thighs or calves.
- + Pain relief: They help to release muscle tension and improve mobility.
- + Fascia training: With regular use, the products can help to prevent muscle soreness and support recovery.

Warranty and liability

The BLACKROLL® fascia products come with a two-year warranty. Damage caused by improper use or failure to follow the care instructions is excluded from the warranty.

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BLACKROLL AG
WSVK Oederan GmbH