MASPO VIBRAMAT

ANLEITUNG

Art.-Nr.  31506 · Kategorie:  Massagegeräte

enjoy your life.

MASPO

Vibramat

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English

Table of Contents

Before use	16
Safety precautions	16
Technical data	16
Basic equipment	17
Optional	17
Standard equipment	17
Special accessories	17
Detailed massage instructions	18
Optimal effects of a complete massage	18
The massage should be in counter-clockwise direction	19
Up to 10 massage variations	19
Start with your right foot	19
Arm and leg massage	20
Unwanted inches	21
Connective tissue massage	21
Neck, Back	21
Muscle massage	22
Sport massage	22
Local massage	22
Internal organs	23
Bronchial asthma and emphysema	23
Lumbago	23
Special care	23
Accessories' upkeep:	23
Chart Massage Schedule	24
Guarantee	25
General Terms and Conditions of Service	26

Before use

No massage in case of feverish colds, varicose veins, thrombosis, phlebitis, pregnancy, jaundice, tumours or other acute ignitions. In case of any kind of illness you always have to consult your physician before starting the massage, when practising special application forms always await the prescription of your physician.

The massage appliance is a wellbeing Product which can help you to feel better. It is not a Product according the medicin-product-law.

Please read this instruction booklet carefully and keep it handy for future reference.

Safety precautions

This appliance conforms to the recognized technical specifications on safety and corresponding safety regulations.

It must only be connected to an AC power supply of the same voltage as that shown in unit identification plate.

Never pull the plug out of its socket by the lead.

Always pull out the mains plug before cleaning or working on the appliance.

This equipment must be disconnected from the mains when not in use.

Children must be supervised during use.

Please do not cover the air intake while using. Overheating-Danger!

Do not by no means plunge the massager into water!

Caution! Do not use this appliance in the bath-tub, shower or above the hand-basin filled with water!

This device has not been designed for being used by persons (including children) that are physically or mentally handicapped or have sensory disabilities. Moreover, it is not allowed to use this device without sufficient experience and/or knowledge, unless these persons will be supervised by persons responsible for their security or in case they have been instructed how to use this device.

Children are to be supervised in order to make sure that they do not play with this device.

Technical data

Motor: 100 W special high-performance motor, interference suppressed, completely insulated, maintenance free. Can be used continuously for several hours.

Revolutions per minute: 1400 - 2800.

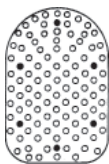
2-speed switch for strong and weak vibration.

Casing and base plate made from shock-resistant plastic material.


Mains lead: CEE-Norm, 2.8 m long.

This unit corresponds to the requirements of 82/499/EC.

Basic equipment


Base plate

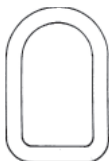Smooth rubber pad

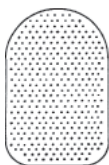ArtNr. 2000 7364

Nylon brush

ArtNr. 2000 7365

Optional

Nodule rubber pad

ArtNr. 2000 7226

Standard equipment

Special accessories

The nylon brush and the smooth rubber pad come with the massager. The nodule rubber pad can be purchased separately.

If you wish to obtain better blood circulation use, according to your sensitivity, one of the two hyperemia brushes.

Hyperamia means: The process by which you improve the blood circulation. The natural bristle brush is more gentle on your skin and is softer than the nylon brush pad.

The smooth rubber pad decreases the massage strokes of the base plate, but at the same time, increases the internal effect of the massage.

To put on the smooth rubber pad or the nodule one, place the straight side of the base plate on the straight side of the pad. Then bend the rounded part of the pad over the rounded front of the bottom plate.

The procedure to put on the nylon brush is the same. First it is important that you adjust firmly onto the ridge the straight end of the brush over the straight part of the base plate.

Then press the rounded part of the brush over the rounded surface, and finally at both sides. In order to take off anyone of the accessories, press your thumbs firmly on both sides of the back of the pad and remove the pad from both hatches.

Detailed massage instructions

It is advisable to spend a few minutes to acquaint yourself with your new large-surface massager. It's to your own advantage since the more you know about the vibratory massage, the more you are going to benefit by it.

The excellent results of the large surface massagers are achieved by their oscillatory effect.

Horizontal vibration plus vertical impulse.

A three-dimensional effect - on the length, the width, and the depth. Muscles and tissue are loosened, the blood circulation is improved and the metabolism is enhanced.

Contrary to a hand massage, the large-surface apparatus does the work 3 to 4 times faster. It even offers another advantage - from the first application, you loose the fear of the often painful hand massage. Large-surface massagers rarely cause pain.

No massage in the kidney area, as this may lead to kidney damages!

The functional and ergonomic design integrates a progressive technologie and offers a maximum efficiency as well as a soft running.

Optimal effects of a complete massage

- **Accessories optional***
- **Duration of massage: 10 -12 minutes**

Do you wish to improve your general health condition?

Then a body or complete massage is suggested once a week.

While a hand massage takes 40 to 50 minutes of your free time, large-surface massager does the work in 10 to 12 minutes.

Do not massage longer!

A lengthy application is useless since the effects obtained could be neutralized by a longer massage.

You can divide the complete massage on two weekdays.

Today, for example, massage your feet, legs, posterior and arms. Three or four days later the body: chest, neck and both side of the spinal cord

* Optional accessories means:

Either without an accessory at all (just the base plate) or with a nylon brush, a smooth rubber pad, a natural bristle brush or a nodule rubber pad.

The massage should be in counter-clockwise direction

Hold the massager by the handle with one hand and with the other hold the knob. Massage slowly and in circles, always in counter-clockwise direction. When massaging the whole body, follow this pattern:

1. Feet
2. Legs
3. Posterior
4. Arms
5. Chest
6. Abdomen
7. Neck
8. Back

If somebody can massage you, lie down and relax your muscles.

The massages should be given with the base plate. Press firm but gently the surface you wish to massage

Strong pressure does not produce a stronger effect, on the contrary, it diminishes the result.

According to your needs and the sensitivity of your skin, you may choose from the four accessories.

Up to 10 massage variations

Do you prefer a stronger or weaker massage effect? You may choose between the two vibratory strengths.

Speed 1 - weak

Speed 2 - strong

And what about the different massage variations?

These can be obtained by using the accessories. You can also use the massager without any accessory, just the base plate.

Start with your right foot

- **Accessory: base plate, nylon brush or smooth rubber pad.**
- **Massage: 1- 2 minutes each foot.**

Start the body massage with your right foot.

At first massage between the two ankles and then from the heels up to the toes. Now you can massage your soles.

Massage your feet only while sitting down.

Do not step on the massager.

An additional massage from the Achilles'tendon to the calf is highly recommended.

Can you feel the results? You may now start massaging your left foot in the same way.

- Do you have to stand while working?
- Do you have swollen feet and legs?

Regular massage helps in many cases. After a few massage strokes you will realize how your tired feet become lighter. The effects of tiredness are vanishing like magic.

It is important that feet are taken care of every day since in most cases the swollen parts slowly disappear.

Flat feet can also be avoided if massaged beforehand. In case of flat feet, the massage will decrease pain, especially combined with active foot gymnastics.

The foot massage with a large-surface massager is gymnastics and massage at the same time.

Do your arms itch?

Do your hands feel numb?

Daily massage can help. Start again with your right hand. If you have heart ailment, massage your left arm vary carefully. Check with doctor first.

- Accessories: base plate.
- Duration: 2 -3 minutes each arm.

Inflexion side

Extensor side

Arm and leg massage

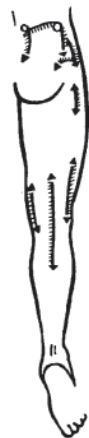
- **Accessories:** base plate, nylon brush pad.
- **Duration:** 2 -3 minutes each leg.

The leg massage improves circulation of the blood and loosens muscular cramps in legs and thighs. Leg cramps can be treated satisfactorily. The calves, the back of your knee and the lower thighs should be also massaged. Even if you have soft skin, this type of massage will be pleasant.

Be careful if you have varicose veins, thrombosis or phlebitis. Ask your doctor beforehand.

Extensor side

Inflexion side

Unwanted inches

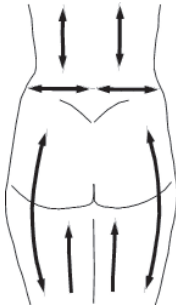
- **Accessories:** see chart massage schedule.
- **Duration:** see chart massage schedule.

Place the massager where you have excessive inches - waist, thighs, hips, posterior, and stomach.

Daily massage will activate the tissue fluid, the retention of water will be reduced, and hence, the loss of "fat" is helped.

Massage as shown below.

For faster results go on a diet simultaneously.

Connective tissue massage

- **Accessories:** optional.
- **Duration:** 2 - 5 minutes.

If you wish to slow down the ageing process you should start on time with the necessary measures to regenerate your cells.

A very positive way to achieve this is by a connective tissue massage. This will bring back the tissue's elasticity.

Metabolic residues slowly disappear as well as the functional problems of the organs involved. Massage as for the complete body massage.

Neck, Back

- **Accessories:** see chart massage schedule.
- **Duration:** 4 minutes approximately.

Long car rides, little movement at work, too much TV watching. We sit down too much. And the results?

Strained neck and shoulder muscles. Back and neck pain.

The large surface massager is excellent for preventive spinal cord massage.

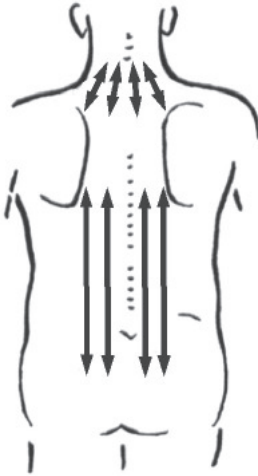
Its oscillatory movement is carried over directly to the spinal cord and cramps in the vertical back muscles will disappear.

Pains here originated can be decreased or disappear.

Massage as shown in the picture.

First the neck (approx. 1 minute). Then right and left of the spinal cord (approx. 1 minute).

i Damaged discs are extremely painful. Preventive vibratory massage with a large-surface massager could avoid this. Check with doctor first.

Muscle massage

- **Accessories:** base plate, smooth rubber pad, noded rubber pad.
- **Duration:** as long as desired.

Cramped or tense muscles?

Due to its three-dimensional massage effect (oscillatory movement) and the large surface of the apparatus you obtain optimum results.

Vibratory massage, especially on the hips, shoulders, and lower back region is very effective to loosen muscular cramps.

Sport massage

- **Accessories:** optional.
- **Duration:** depending upon which part of the body, 1- 4 minutes.

Not only the muscles directly strained by sport are massaged, but also those stimulated indirectly.

The large-surface massager is also recommended for a short massage during a championship of training since not only is the massage intensive, but most important of all, it can be performed during a short period of time.

A decisive advantage:

You may massage without removing your clothes.

Local massage

- **Accessory:** smooth rubber pad.
- **Duration:** 1- 2 minutes.

Vibratory massage not only relaxes and keeps you fit but also can be employed satisfactorily for local massage.

The period of malaise can very often be shortened. Its very differential oscillating movement obtains even better results than a massage by hand.

Warning:

Use carefully. Consult a physician before you use it for local massage.

Further information can also be obtained at book shops.

Internal organs

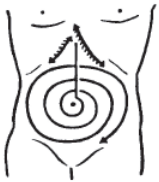
- **Accessory:** smooth rubber pad.
- **Duration:** 1- 2 minutes.

Organs can be positively influenced and the metabolism improved by massaging a specific part of your body.

Ask your doctor, masseur or gym therapist for these special cases.

To stimulate bowel movement, massage your back and stomach as shown on the two pictures below. When massaging your stomach, use the slow speed.

No massage in the kidney area, as this may lead to kidney damages!

Front

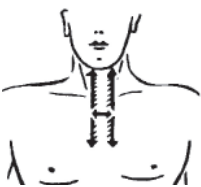
Back

Bronchial asthma and emphysema

- **Accessories:** smooth rubber pad.
- **Duration:** according to doctor's advice.

Relief is obtained by massaging slowly and smoothly. Massage both sides of the larynx, the breastbone, and also between your spinal cord and shoulders. Check with your doctor first.

Front

Back

Lumbago

- **Accessories:** base plate, noded rubber pad.
- **Duration:** 3 - 5 minutes.

This sudden back and neck pain needs a daily massage.

Massage the lower back at both sides of the spinal cord.

Check with your doctor first.

Special care

This large-surface massager is absolutely maintenance free.

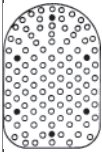
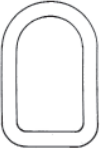
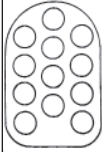
Accessories' upkeep:

The accessories should be cleaned from time to time by washing them in lukewarm soapy water.

The base plate which is permanently fastened to the massager can be cleaned with a damp cloth. DO NOT immerse the massager in water under any circumstances.

Chart Massage Schedule

Duration of single-massage,
vibration speed and
instructions for use

Parts of the body	Duration in minutes	Speed				
arms	2 - 3	I II	●			●
stomach	1 - 2	I			●	
legs	2 - 3	II	●	●		●
connective tissue	2 - 5	I II	●	●	●	●
chest	1 - 2	I			●	
double chin	up 1	I			●	
between calves and feet	1 - 2	II			●	
feet	1 - 2	II	●	●	●	
body*	10 - 12	I II	●	●	●	●
posterior	2 - 3	II	●	●		●
face	up 1	I			●	
hands	1 - 2	I II	●			●
lumbago	3 - 5	I II	●			
hips	2 - 3	II	●	●	●	
organs	1 - 2	I			●	
lumbar region	1 - 2	I			●	
muscles	1 - 5	I II	●		●	
neck	1	I II	●		●	
thighs	2 - 3	II	●	●	●	●
back	2	II				●
shoulders	2	II	●			
waist	1 - 2	II		●	●	●
calves	2 - 3	II	●	●		
spinal cord relax	1	I II			●	

*does not equal individual massage duration for each part of the body a full massage can be performed quicker.

Please keep this address in a safe place together with the installation guide.

To help us answer your questions quickly and competently please provide the information printed on the type shield including the model, item no. and serial no., in all inquiries.

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www.eos-sauna.de

WARRANTY

The warranty is provided according to the legal regulations at present.

Manufacturer's guarantee:

- The period of guarantee starts from the date of purchase and lasts up to 2 years by commercial use and 3 years by private use.
- Always include the completed guarantee certificate when returning equipment.
- The guarantee is void for appliances which have been modified without manufacturer's explicit agreement.
- Damages caused by incorrect operation or handling through non-authorized persons are not covered under the terms of guarantee.
- In the event of a claim please indicate the serial number as well as the item number and model name with detailed description of the fault.
- This guarantee covers defective parts and labour but not the defects caused by wear and tear.

In case of complaint please return the equipment in its original packaging or other suitable packaging (caution: danger of transport damage) to our service department.

Always include the completed warranty certificate when returning equipment.

Possible shipping costs arising from the transport to and from point of repair cannot be overtaken by us.

Outside of Germany please contact your specialist dealer in case of warranty claims. Direct warranty processing with our service department is in this case not possible.

Equipment commissioning date:

Stamp and signature of the authorized electrician:

General Terms and Conditions of Service

I. Scope

Unless otherwise agreed in writing in a specific case, these terms and conditions of service shall apply to service operations, including examining and repairing complaints. All our existing or future legal relationships shall be governed solely by the following terms and conditions of service. Our recognition of any conflicting terms and conditions of the Ordering Party shall be conditional upon our having given our express written consent to their applicability. We hereby expressly object to any terms and conditions of the Ordering Party contained in its General Terms and Conditions of Business or order confirmation. If order confirmations or deliveries are accepted without reservation, this shall not be deemed to constitute recognition of such terms and conditions. Any ancillary agreements or amendments must be confirmed in writing.

II. Prices and Payment Terms

The Ordering Party shall bear the following costs in connection with the service operation:

- Disassembly/assembly and electrical (de-) installation
- Transportation, postage and packaging
- Function testing and troubleshooting including inspection and repair costs

There shall be no third-party billing.

III. Performance Obligations / Ordering Party's Cooperation

The Ordering Party shall provide free-of-charge assistance to the manufacturer in carrying out the service operation.

In the case of a warranty claim, the manufacturer shall make the replacement parts necessary for the service available to the Ordering Party free of charge.

IV. Service Visit by the manufacturer

In the event that it is essential that a manufacturer employee carry out the service operation on site, this must be agreed in advance. Where the main reason for the service call is not the fault of the manufacturer, any costs incurred shall be recharged to the Ordering Party after the service visit.

V. Liability

The manufacturer shall assume liability in accordance with the currently applicable statutory regulations. The packaging for all of our products is designed for the shipping of individually packed goods (pallet). We expressly point out that our packaging is not suitable for individual

shipments via parcel post. The manufacturer shall accept no liability for damage incurred as a result of improper packaging in an individual shipment.

VI. Manufacturer's Warranty

The manufacturer's warranty shall apply only in the event that installation, operation and maintenance have been carried out in accordance with the manufacturer's specifications contained in the assembly instructions and instructions for use.

- The warranty period shall commence from the date on which proof of purchase is provided and shall be limited, in principle, to 24 months.
- Warranty services shall be performed only if the proof of purchase relating to the equipment can be presented.
- Any and all warranty claims shall become void if modifications are made to the equipment without the manufacturer's express consent.
- Any warranty claim shall likewise become void in the case of defects that arise due to repairs or interventions made by unauthorized persons or due to improper use.
- In the case of warranty claims, the serial and article numbers must be indicated together with the designation and a meaningful description of the fault.
- This warranty shall cover defective equipment parts, with the exception of normal wear parts. Wear parts shall include, among other things, light sources, starters, gas or oil pressure dampers as well as acrylic glass panes, tubular heating elements and sauna heater stones.
- Only original replacement parts may be used within the warranty.
- Service visits by outside companies shall require a written order to be issued by our service department.
- The equipment in question shall be sent to our service department by the Ordering Party and at its expense.
- Electrical assembly and installation work, including in the event of service or replacement, shall be carried out at the Customer's expense and shall not be borne by the manufacturer.

Complaints in respect of our products shall be reported to the specialist trader responsible and shall be exclusively handled via the latter.

The manufacturer's General Terms and Conditions of Business, as amended, shall apply in addition to the foregoing terms and conditions of service.

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Sport-Tec GmbH