



GLEICHGEWICHTSSTÖCKE FÜR SKITRAINER PRO FITTER 2 ANLEITUNG



3050-2600 Portland St. SE, Calgary, Alberta, Canada T2G 4M6
Corporate Office: 1-800-348-8371 fax: 1-866-250-8824 e: sales@fitter1.com
www.fitter1.com • 1-800-fitter-1

BALANCE AID PRODUCT ASSEMBLY & WARRANTY

Fitter International Inc. guarantees your Balance Aids against manufacturing defects. If your Balance Aids are received in a defective state, please contact us at **1-800-fitter-1** or **sales@fitter1.com** a customer Sales Representative will be happy to assist you.



Your Box Should Contain:


- 2 Aluminum Poles
 - Rubber mounts attached
 - Screws attached to mounts
- 2 Touch Cap Spheres
- 2 Large 8 inch Flat Discs
- 1 Allen Key

1



Using the Allen Key remove screw from rubber mount. Fit screw through hole in centre of large 8 inch Flat Disc. Insert the exposed end of the screw back into the rubber mount.

2

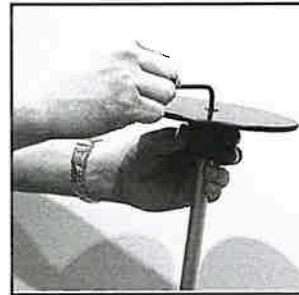

Again using the Allen Key, tighten the screw to lock the flat disc in place. Make sure the screw is snug, but do not over-tighten.

3



Place the Touch Cap Sphere snugly onto the top of the pole. **Please note the Touch Cap Sphere is dual purpose and is force fit onto the aluminum pole, there is no nut on the top.**

4



Very Important: Firmly hit the sphere with a hammer or turn upside down and hit firmly against the floor to lock into place.

Thank you for choosing Balance Aids as part of your rehabilitation or fitness program. We appreciate the opportunity to serve you and hope to hear from you again soon. Visit **www.fitter1.com** for our complete line of quality fitness products.

Vertrieb durch:



==== www.sport-tec.de ====

Sport-Tec GmbH
Physio & Fitness
Lemberger Str. 255
D-66955 Pirmasens

Tel.: +49 (0) 6331 1480-0
Fax: +49 (0) 6331 1480-220
E-Mail: info@sport-tec.de
Web: www.sport-tec.de